


READ ME FIRST

Installing Spartan 06 for Mac OS X

To install Spartan 06 on your hard drive, insert the Spartan 06 CD into your CD-ROM drive. Open the Spartan 06 CD by double clicking on its icon, and drag the "Spartan 06" icon to your hard disk. (Dragging "Spartan 06" to your "Applications" folder is recommended.)

Installing the Spartan Database

Spartan 06 utilizes a database of common molecular structures (on a separate DVD named "Spartan Database"). This data (approximately 3.1 GB) consists of three folders:

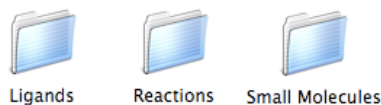

To make your database, create a "Spartan Database" folder on your hard disk, preferably in the "Applications" folder. Then drag the three folders from the "Spartan Database" DVD to your new database folder. Spartan 06 will ask you to locate these folders when it starts for the first time (below).

There are two other data folders on the Spartan Installation CD:



The "Library" folder contains collections of useful molecular fragments that can be added to the "Custom" pop-up menu in the Model Kit (see p. 200 in the User's Guide). It is suggested that you drag this folder into your "Spartan Database" folder as well.


Finally, all of the completed tutorial exercises from the User's Guide can be found in the "Tutorials" folder.

Spartan 06 Licensing

Spartan 06 supports three licensing methods:

- Local USB hardware keys (a small device attached to the computer's USB port)
- NetHASP licensing (the computer "checks out" a license from a pool available on the network)
- Access codes (which are time-limited for demonstration/evaluation purposes).

If you are using a local USB key, you will need to install drivers to support the key before starting Spartan 06 for the first time. **UNPLUG ANY LOCAL HASP KEYS** attached into this computer, then locate the "Hardware Key Drivers" folder on the CD and double click the "HASP Universal Installer" icon. A window will be displayed:

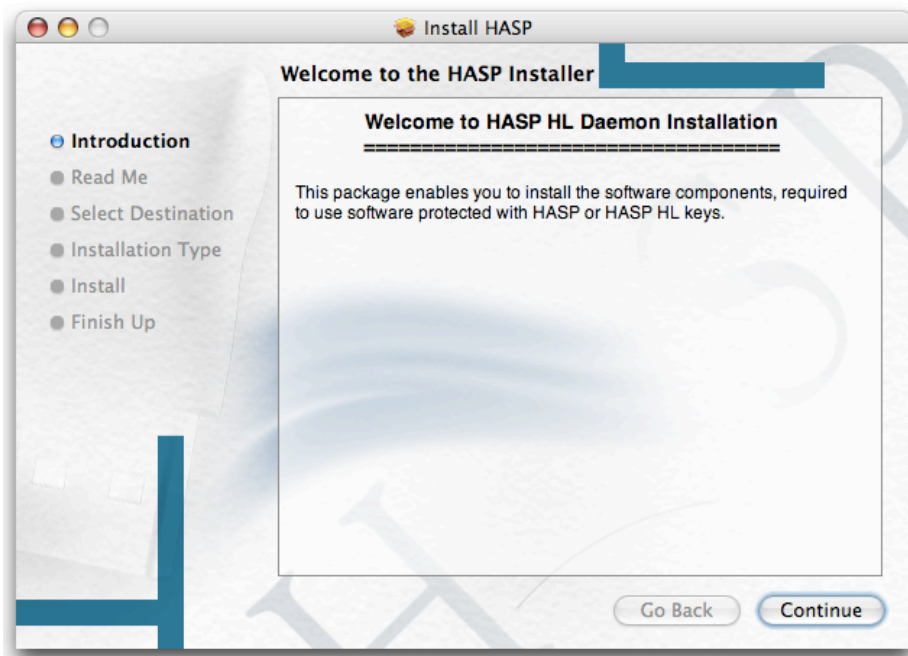


Click Continue, entering your password when prompted. Quit the Installer when installation is complete, re-insert any local hardware keys, and you will be ready to start Spartan '06.

If you are using a demonstration access code, you will need to connect the computer to the Internet when you enter the code to activate the temporary license. (After activation, Spartan doesn't need an internet connection unless you want to perform Protein Data Bank searches- see page 314 of the manual.)


Launching and Licensing Spartan '06

When Spartan '06 is started for the first time, it will display an setup window:

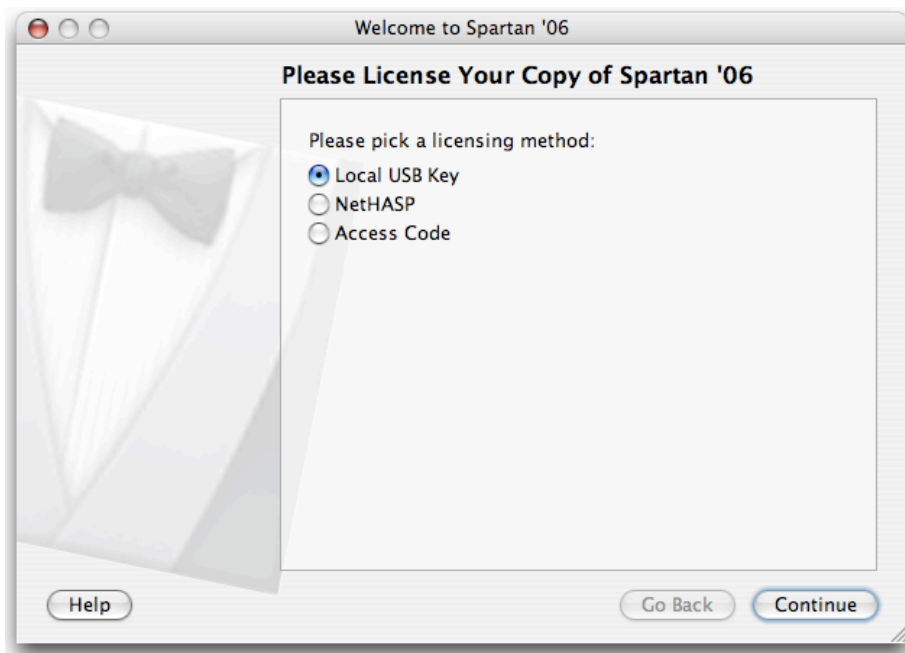

Select your licensing method. If you are using a demonstration access code, another window will be displayed to enter the access code sequence and validate the code over the internet.

After Spartan has been licensed successfully, you will be prompted to locate the folders that you dragged to your Database folder:

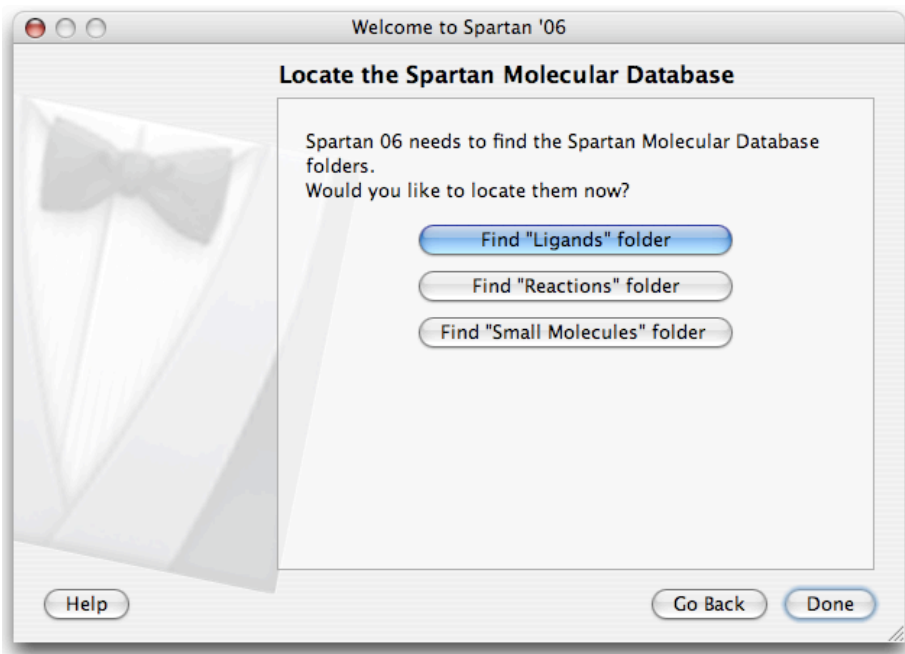

Click the buttons and locate the folders in the file browser that appears. When you have locate all three folders, click "Done" and you will be presented with the Spartan '06 main window.